



# 1100 - Wireless Temperature Sensor

Revision: 3.0 | DS155

## Table of contents


---

Introduction .....	3
Pack Contents .....	4
Operational Overview .....	5
Connectivity .....	8
Charging the Sensor .....	9
Firmware Updates .....	10
Usage Information .....	11
Practical Investigations .....	13
Sensor Specifications .....	15
Limited Warranty .....	16
Compliance .....	17
Troubleshooting .....	18
Notices .....	19
Contact Information .....	20
PDF Translations .....	21

## Introduction

---

Thank you for purchasing the Smart Wireless Temperature Sensor. We pride ourselves on producing high quality products that meet with the demands of the busy classroom environment. If you have any problems using this sensor, please read this documentation in full before contacting the Data Harvest support team.

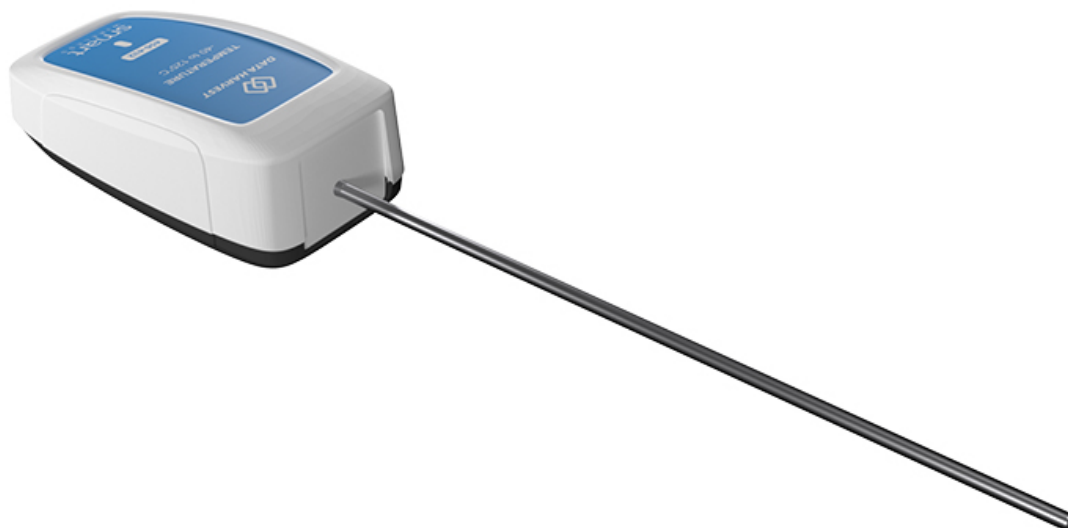

## Overview

The Smart Wireless Temperature sensor has a thermistor housed at the end of a stainless steel tube. The tube is 3 mm by 160 mm and is made from AISI 316 stainless steel, which has good thermal transfer and high corrosion resistance. The sensor is suitable for use with experiments in chemistry, biology, physics, earth & environmental science.

The Smart Wireless Temperature sensors are USB and Bluetooth compatible. Using Bluetooth, a sensor can connect to mobile devices, tablets laptops and desktops.

The sensors are supplied with a mini USB lead (1 m standard A to standard mini B) and calibrated with degrees Celsius ( C) as the default range.

---

## Pack Contents

---

This product is supplied with the following items:

- [1 x Wireless Temperature Sensor](#)
- 1 x USB Connecting Lead


## Additional Accessories

To get the most from your Smart Wireless Fast Response Temperature Sensor, the following items should be considered:

- [Wireless Fast Response Temperature Sensor](#)
  - [Wireless Infrared Sensor](#)
-

## Operational Overview

The diagram below shows the specific parts of the sensor. Read further to explore the functionality of each part of the sensor.

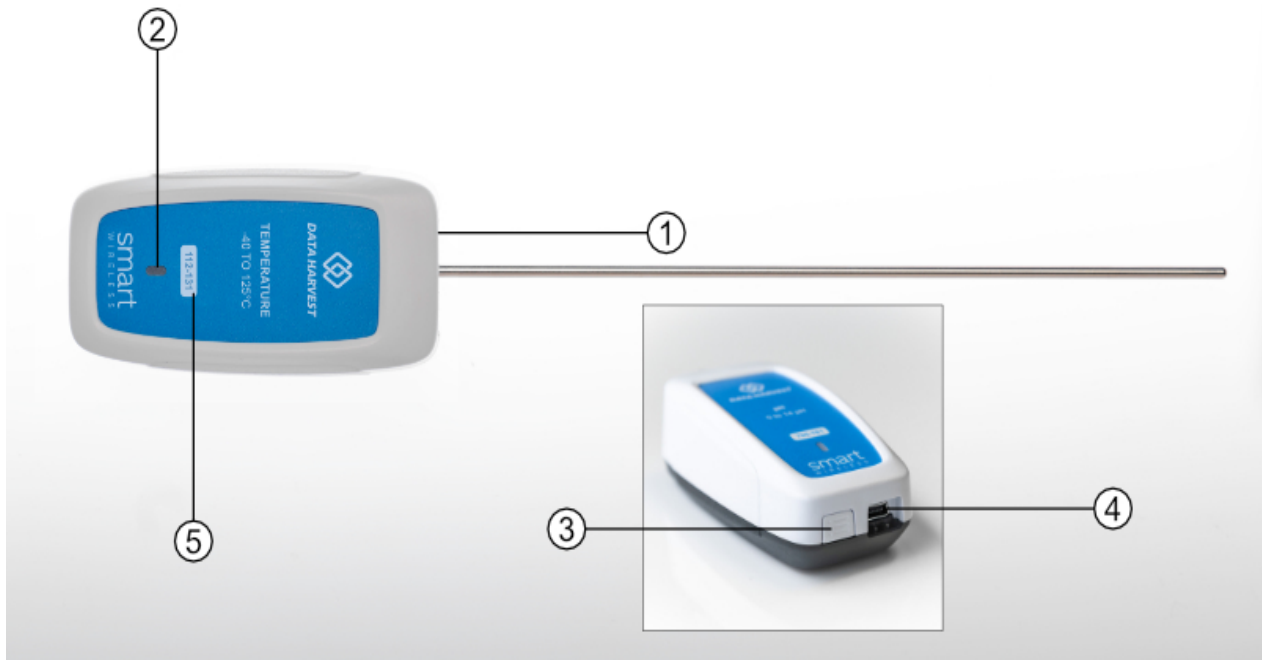

1. Sensor End Cap
2. Status Indicator
3. On/Off Switch
4. USB Port
5. Unique ID Number







### Sensor End Cap (1)

Most Smart Wireless Sensors feature an end cap that is specific to the requirements of the device's internal sensor. The sensor's end cap is the direct interface between the device's internal sensor and your experiment.

### The Status Indicators (2)

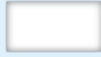The sensor features a single status indicator that changes colour and flashes. See the table below for further information.

Status Light	Indicates
No light	Sensor is Off. Short press the On/Off switch
Blue flashing	Sensor is On and Bluetooth advertising

White flashing		Charging via USB mains charger or USB port, Sensor is On and Bluetooth advertising
Red, Green, Blue Flashing		Charging via USB mains charger or USB port, Sensor is Off
Green flashing		Communication with the EasySense app (via USB or Bluetooth) has been established
Solid Green		Fully charged
Orange flashing		Recording data, a fast pulse indicates awaiting trigger in Remote Mode
Red flashing		Battery is low

### On/Off Switch (3)

The sensor's on/off switch allows you to turn the sensor on, off or perform a hard reset.

#### To switch the sensor off

- Press and hold down the On/Off switch until the white light shows, then release.
- If not communicating with the EasySense app, the sensor will turn off after a period of one hour of inactivity.

#### Hard resetting the sensor

- If necessary, attach the sensor to power.
- Press and hold down the On/Off button for at least 8 seconds until the status LED gives a flash of blue light, then release.
- If the sensor fails to respond, contact Product Support at Data Harvest. Please provide details of:
  - The computer platform it is being used with and the EasySense app's version number.
  - A description of the problem being encountered.

### USB Port (4)

Use to connect to a computer or a charging unit.

For specific USB or Bluetooth connectivity instructions, please see the 'Connectivity' section of this documentation.

For instructions on charging your device, see the section on 'Charging the Sensor'.

### Unique ID Number (5)

All Smart Wireless Sensors are labelled with a unique ID number. This number is used in the EasySense app, so that you can identify each sensor when making a connection wirelessly.

### The Sensor and EasySense

Please make sure that you use the latest release of the EasySense series of software. Both collection and analysis of data is available here, on a variety of operating systems.

#### Direct Data Logging

The sensor is designed to work directly with EasySense (as an installed application or PWA). A full compliment of experiments can be run by using the sensor through Bluetooth™ or USB. EasySense will support direct logging and data storage when connected as above.

**Remote Data Logging**

The ability to capture data independently (free of a capture station) is done through EasySense's Remote Mode.

This facility may be found in EasySense, under Setup. Once the conditions for data collection have been established, the sensor can be set to initiate collection for example, using a rapid press of the power button. Initiation of the experimental data collection by the software is followed by remote detachment; collection is then on the sensor.

Data gathering is realised by using Setup once again.

Details are given in the latest EasySense User Guide.

---

## Connectivity

---

The sensor is both USB and Bluetooth compatible. Install the EasySense app, if it is not already on your device. For details of how to operate the EasySense app, please refer to the EasySense documentation.

### USB Connectivity

#### Quick Steps

1. Connect the sensor to the computer's USB port using the USB cable supplied.
2. The computer will automatically detect a new device and depending on your operating system, will install any applicable device drivers.
3. Start EasySense app.
4. Within the EasySense app, the Devices icon will change to green to show that the sensor is connected, and the status light on the sensor will also turn green.
5. Begin your practical investigations.

### Bluetooth Connectivity

Using Bluetooth, the sensor can wirelessly connect to mobile devices such tablets and mobile phones, as well as desktop or laptop computers, giving students the ability to run experiments independently without being tethered to a device.

See the EasySense app user manual system requirements for further details.

#### Quick Notes on Bluetooth Connectivity

Only use with the EasySense app, you do not need to pair the device. If paired, the sensor will not be available to the EasySense app.

Computers or devices will need to support Bluetooth Low Energy (BLE). For further information refer to the instructions provided for the EasySense app.

#### Quick Steps

1. Short press the on/off switch to turn the sensor on, blue LED will flash.
  2. Open the EasySense app.
  3. Select the Devices icon.
  4. Select your sensor from the list of available sensors to connect to the device. Your sensor is identified by its unique ID in the list.
  5. Click on connect at the side of your sensor in the list.
  6. The Devices icon will change to green and the status light on the sensor will flash green to indicate a connection has been established.
  7. Begin your practical investigations.
-



## Charging the Sensor

---

The Smart Wireless sensors are fitted with a rechargeable lithium-ion battery and can be charged via the USB port. Use the supplied USB lead to connect the sensor either directly to a USB port on your computer, a powered USB hub or a USB mains charger that outputs 5 V at 500 mA or more.

A full charge can take up to four hours.

### Additional Information

Whenever the sensor is connected to the USB port on the computer or to a USB mains charger (output 5 V at 500 mA or more), it will automatically recharge the battery (LED status flashing white).

When connected to a computer, the computer should be turned on and not in sleep or standby mode, as the battery may drain instead of charge.

The sensor will stay awake for five minutes when Bluetooth advertising (LED status flashing blue).

Lithium-ion batteries are 'memory-free' and prefer a partial rather than a full discharge. Constant partial discharges with frequent recharges will not cause any harm. Frequent full discharges should be avoided whenever possible. Ideally the sensor should be stored at about 40% or more charge.

The speed at which a lithium-ion battery will age is governed by both its storage temperature (preferably less than 40 C) and state-of-charge.

---

## Firmware Updates

---

Occasionally Data Harvest may release updated firmware which will contain improvements or new features.

Updates will take place when you connect your sensor to the EasySense app. You will be given the option to decline an update.

Updates can be performed over USB or Bluetooth and will typically take less than one minute. Updating firmware over USB will be quicker than Bluetooth.

Do not disconnect the sensor, or power off during the update.

If you have a wireless connection to the EasySense app, the sensor will have to be reconnected after performing the update.

---

## Usage Information

---

The thermistor sensor is at the tip of the stainless steel tube. The thermal transfer to the tube, and then to the thermistor, could make surface temperature measurement less accurate. If you need to measure the surface temperature, place the tube along the surface; probing with the tip will not allow enough heat to transfer.

Temperatures higher than 150°C (302°F) will damage the sensor.

Do not put any part of the sensor into a direct flame or onto a hot plate.

The electronics in the Smart Wireless module of the sensor have a working range of 0 to 40°C, so cannot be subjected to the same extremes of temperature as the stainless steel portion. Position the sensor so that this part of the sensor is kept away from the source of heat. Take steps to protect against excessive steam.

The Smart Wireless module is not waterproof; it has an operating range of 0 to 95% RH (non-condensing). Do not place the sensor in an environment in which high humidity levels are possible, as this may result in damage or malfunction.

Do not submerge the sensor beyond the stainless steel portion.

The metal tube section can be washed after use. Clean the Smart Wireless module using a damp cloth; do not immerse in water or detergent. If the sensor has been left in the cold, let it warm to near room temperature before waking it from sleep.

AISI 316 stainless steel has a high resistance to corrosion from a wide variety of weak acid and alkalis. Some environments e.g. saltwater, may cause some discolouration to the stainless steel tube but this will have no effect on the sensor's performance.

The sensor can be left in an alkaline solution such as sodium hydroxide (NaOH) for up to 48 hours, with only minor discolouration. We do not recommend use in a solution whose concentration is greater than 3 mol dm<sup>-3</sup>.

The maximum length of time recommended for exposure to an acid is dependent on the acid's concentration. In general, we do not recommend that Temperature sensors be left to soak in acids of between 1 to 3 mol dm<sup>-3</sup> concentrations for longer than 48 hours.

The exceptions to this rule are hydrochloric acid (HCl), and sulphuric acid (H<sub>2</sub>SO<sub>4</sub>).

### The maximum recommended times for exposure are:

HCl	Maximum Exposure
1 mol dm <sup>-3</sup>	20 Minutes
2 mol dm <sup>-3</sup>	10 Minutes
3 mol dm <sup>-3</sup>	3 Minutes

H <sub>2</sub> SO <sub>4</sub>	Maximum Exposure
1 mol dm <sup>-3</sup>	48 Hours
2 mol dm <sup>-3</sup>	20 Minutes
3 mol dm <sup>-3</sup>	10 Minutes



## Practical Investigations

The Smart Wireless Temperature Sensor can be used to investigate a number of scientific experiments such as:

- Monitoring indoor and outdoor temperature
- Weather studies
- Insulation studies
- Solar homes
- Monitoring endothermic and exothermic reactions
- Solubility of salts
- Studying freezing and boiling points
- Cooling rates
- Evaporation, radiation, conduction and convection investigations
- Energy content of fuels and foods
- Heat of fusion investigations

## Online Videos

Learn how to use data logging in the classroom with our Secondary Science Academy demonstration videos, which will walk you through using the new EasySense app and show you how to get hands-on with the latest Bluetooth wireless sensors. The video experiments will show you how to get the best out of your science lessons.

New online content is being continuously uploaded onto our YouTube channel, including practical worksheets as well as videos.


See our website for further information and links.



### Explore Bluetooth Sensors

Are you looking to make the jump to our smart wireless sensors? Or have you recently purchased them and want to know more about how they work?

[View video playlist](#)

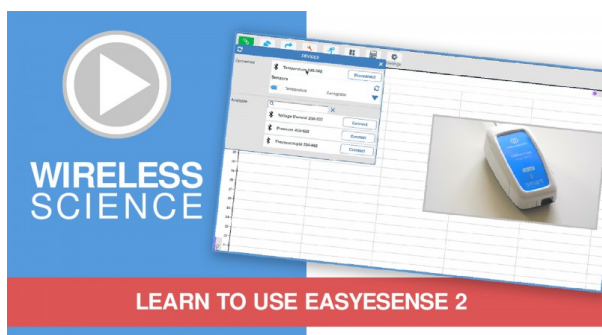

### Explore EasySense

The core of our science platform is our EasySense app. In these videos you will learn everything from the basics of our software to the most in-depth features.

[View video playlist](#)



## Explore Science Practicals

See our Smart Wireless Sensors in action with a range of practical experiments. This is the best way to get started with the new Bluetooth sensors!

[View video playlist](#)

## Sensor Specifications

---

Please read the following table for sensor specifications.

Feature	Detail
Measurement Ranges	-40 C to 125 C and Fahrenheit: -40 F to 257 F
Accuracy	±0.3 C (0.5 F) at 0 to 70 C, rising to ±0.6 C (1.1 F) at extremes of range
Resolution	0.01 C (0.01 F)
Fastest logging speed	50 samples per second [20 ms]
Connectivity	Wired via USB Wireless via Bluetooth
Bluetooth Specifications	Bluetooth 4.2 low energy radio, single mode compliant Transmit (TX) power: 0 dBm Receiver (RX) sensitivity: - 90 dBm Usable transmission range: up to 10 m in open air Frequency Range: 2.402 to 2.480 GHz operation Operating range: 0 to 40 C and 0 to 95% RH (non-condensing)
Internal Battery	Rechargeable internal lithium-ion 3.7 V
Storage/Operating Temperature	0 to 40 C
Humidity	0 to 95% RH (non-condensing)
Physical Specifications	Weight: approx. 80 g External dimensions: approx. height 33 mm x width 50 mm x length 90 mm Stainless steel rod: 3 mm dia. x 160 mm

---

## Limited Warranty

---

For information about the terms of the product warranty, see the Data Harvest website at: <https://data-harvest.co.uk/warranty>

## Product Repairs

When returning goods to Data Harvest, please download and complete the repair return [form](#) to ensure you have sent us all the information we require, and send it to us alongside the item to be repaired. The second page of this form includes a return address label.

If you have purchased a Data Harvest manufactured product via a different company, please also supply proof of purchase.

## Postage Charges

- In the event of a fault developing, the product must be returned in suitable packaging to Data Harvest for repair or replacement at no expense to the user other than postal charges.
- There will be no postal charge for the return of repaired goods to any mainland UK address (for other areas, additional shipping charges may apply).

## Out of Warranty Repairs

Please visit <https://data-harvest.co.uk/repairs> for the most up to date charges for out of warranty repairs.

## Warranty on Repaired Items

Once an item has been serviced and repaired, the product will have 1 year warranty against further failure of the component repaired.

## International Returns

Please contact the authorised Data Harvest representative in your country for assistance in returning equipment for repair.

---



## Compliance

---

This product complies to the following standards:

### Waste Electrical and Electronic Equipment Legislation

Data Harvest Group Ltd is fully compliant with WEEE legislation and is pleased to provide a disposal service for any of our products when their life expires. Simply return them to us clearly identified as 'life expired' and we will dispose of them for you.

### FCC Details


This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

### CE

This product conforms to the CE specification. It has been assessed and deemed to meet EU safety, health and environmental protection requirements as required for products manufactured anywhere in the world that are then marketed within the EU.

### UKCA

This product conforms to the UKCA specifications.

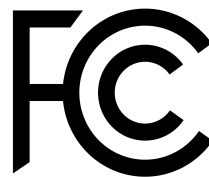

## Troubleshooting

---

If you experience any problems with your product, please try the following troubleshooting tips before contacting the Data Harvest support team.

Feature	Detail
Loss of Bluetooth Connectivity	<p><b>If the sensor loses Bluetooth connection and will not reconnect try:</b></p> <p>Closing and reopening the EasySense app.</p> <p>Switching the sensor Off and then On again.</p> <p>If you are using a Bluetooth Smart USB Adaptor on your computer, unplug the adaptor, plug back in again and try to reconnect.</p> <p>Hard reset the sensor and then try to reconnect.</p>

---

## Notices

---

Please read the following notices with regards to using your sensor

1. The sensor is much smarter than traditional Bluetooth sensors and you are not required to pair the device. If paired, the sensor will not be available to the EasySense app.
  2. When the sensor is connected to a computer, the computer should be turned on and not in sleep or standby mode as the battery may drain instead of charge.
  3. Data Harvest products are designed for educational use and are not intended for use in industrial, medical or commercial applications.
  4. The sensor is not waterproof.
  5. Plastic parts may fade or discolour over time if exposed to UV light. This is normal and will not affect the operation of the sensor.
-

## Contact Information

---

To contact Data Harvest directly, please use any of the following channels:

### Traditional Communications

Data Harvest Group Ltd.  
1 Eden Court, Eden Way,  
Leighton Buzzard,  
Bedfordshire,  
LU7 4FY  
United Kingdom

**Tel:** +44 (0) 1525 373666

**Fax:** +44 (0) 1525 851638

**Sales email:** [sales@data-harvest.co.uk](mailto:sales@data-harvest.co.uk)

**Support email:** [support@data-harvest.co.uk](mailto:support@data-harvest.co.uk)

### Online Communications

We have active social media support channels using the following platforms

- [Facebook](#)
- [X](#)
- [YouTube](#)

### Office Opening Hours

Monday to Thursday - 08:30 to 16:45

Friday - 08:30 to 13:30

Saturday & Sunday & UK Bank Holidays - Closed

---

## PDF Translations

---

The PDF formatted download of this manual is by default provided in the English (United Kingdom) language. If an alternative translation is available, it will be listed here.

We have for your convenience included a webpage translation feature to the online documentation which will allow you to translate and print individual pages of this documentation.

---